



IN THE UNITED STATES PATENT AND TRADEMARK OFFICE

Application of: Vincent PETIARD *et al.*

Confirmation No.: 1875

Application No.: 09/849,139

Group Art Unit: 1637

Filing Date: May 4, 2001

Examiner: Joyce Tung

For: METHOD OF DETERMINING THE GENETIC  
MATERIAL OF COCOA IN FERMENTED OR  
ROASTED BEANS AND CHOCOLATE

Attorney Docket No.: 88265-4022

**SUBMISSION OF PRIORITY DOCUMENT**

**Mail Stop Amendment**

Commissioner for Patents

P.O. Box 1450

Alexandria, Virginia 22313-1450

Sir:

Applicant has claimed priority under 35 U.S.C. § 119 to European Patent Application No. 98121043.8, filed on November 5, 1998. In support of this claim, a certified copy of this application is submitted herewith.

No fee or certification is believed to be due for this submission. Should any fee be required, however, please charge such fee to Winston & Strawn LLP Deposit Account No. 50-1814.

Respectfully submitted,

Date 9/1/04

Rodney J. Fuller  
Rodney J. Fuller (Reg. No. 46,714)  
For: Allan A. Fanucci (Reg. No. 30,256)

**WINSTON & STRAWN LLP**

Customer No. 28765

202-371-5904

**BEST AVAILABLE COPY**

THIS PAGE BLANK (USPTO)



**Europäisches  
Patentamt**

**European  
Patent Office**

**Office européen  
des brevets**

**Bescheinigung**

**Certificate**

**Attestation**

Die angehefteten Unterlagen stimmen mit der ursprünglich eingereichten Fassung der auf dem nächsten Blatt bezeichneten europäischen Patentanmeldung überein.

The attached documents are exact copies of the European patent application described on the following page, as originally filed.

Les documents fixés à cette attestation sont conformes à la version initialement déposée de la demande de brevet européen spécifiée à la page suivante.

**Patentanmeldung Nr.    Patent application No.    Demande de brevet n°**

98121043.8

**CERTIFIED COPY OF  
PRIORITY DOCUMENT**

Der Präsident des Europäischen Patentamts;  
Im Auftrag

For the President of the European Patent Office

Le Président de l'Office européen des brevets  
p.o.

**R C van Dijk**

DEN HAAG, DEN  
THE HAGUE,    21/06/04  
LA HAYE, LE

**THIS PAGE BLANK (USPTO)**



Europäisches  
Patentamt

European  
Patent Office

Office européen  
des brevets

**Blatt 2 der Bescheinigung**  
**Sheet 2 of the certificate**  
**Page 2 de l'attestation**

Anmeldung Nr.:  
Application no.:  
Demande n°: **98121043.8**

Anmeldetag:  
Date of filing: **05/11/98**  
Date de dépôt:

Anmelder:  
Applicant(s):  
Demandeur(s):  
**SOCIETE DES PRODUITS NESTLE S.A.**  
**1800 Vevey**  
**SWITZERLAND**

Bezeichnung der Erfindung:  
Title of the invention:  
Titre de l'invention:

**Use of DNA identification techniques for the determination of genetic material of cocoa in  
fermented or roasted beans and chocolate**

In Anspruch genommene Priorität(en) / Priority(ies) claimed / Priorité(s) revendiquée(s)

Staat:  
State:  
Pays:

Tag:  
Date:  
Date:

Aktenzeichen:  
File no.  
Numéro de dépôt:

Internationale Patentklassifikation:  
International Patent classification:  
Classification internationale des brevets:

**C12Q1/68, // A23G1/00**

Am Anmeldetag benannte Vertragsstaaten:  
Contracting states designated at date of filing:  
Etats contractants désignés lors du dépôt:

**AT/BE/CH/CY/DE/DK/ES/FI/FR/GB/GR/IE/IT/LI/LU/MC/NL/PT/SE**

Bemerkungen:  
Remarks:  
Remarques:

THIS PAGE BLANK (USPTO)

5 Nov. 1998

**Use of DNA identification techniques  
for the determination of genetic material of cocoa  
in fermented or roasted beans and chocolate**

The present invention relates to the use of DNA identification techniques for the determination of genetic material of cocoa in fermented and/or roasted beans and chocolate. In particular, the present invention pertains to the use of techniques selected from the group consisting of PCR, RAPD, RFLP or microsatellite identification for the determination of cocoa varieties and/or cacao varieties that have been modified by common breeding techniques or that have been modified by genetic engineering in processed beans and chocolate.

Due to the globalisation of the world a huge variety of food products from all over the world have become available in local marketplaces. This also applies to the raw materials or intermediate products used in the various stages of food production.

Apart from new foreign food products entering local markets, novel, modified products may be obtained, that show a variety of different traits. To this end, e.g. tomato paste deriving from genetically engineered tomatoes have been sold in the United Kingdom since 1996 after having been imported from the United States of America.

Due to national regulations some food products and raw materials therefore have to be labelled so as to be able to determine their origin and/or quality level. One of the drawbacks of such labelling procedures resides in that said labels may accidentally or intentionally be exchanged so that the purchaser may not always fully rely on this sort of identification.

Therefore, there exists a need in the art for a reliable method so as to determine particular raw materials in food products of interest or their origin.

With the progress of DNA analysis techniques efficient methods for the identification of the genetic origin of raw materials/food products have become accessible. However, these techniques may only be applied to products/raw materials the genetic material of which has not been degraded to a substantial degree.

In this respect, R. Greiner et al. report in "Is there any possibility of detecting the use of genetic engineering in processed foods", *Z. Ernährungswissenschaft* **36** (1977), 155-160, of the option that some food products, such as pizza tomatoes, fried potatoes etc. allow for the detection of their genetic material by means of PCR (Polymerase Chain Reaction; Ehrlich et al., *Science* **253** (1991), 1643-1651). In the same paper (Greiner et al., *supra*) it is, however, also mentioned that some processed food products, such as tomato soup, tomato starch, mashed potatoes etc. do not enable such an identification, due to the DNA thereof obviously being degraded to an extensive degree.

For cocoa it was generally accepted that the process steps involved in the preparation of cocoa, such as bean fermentation, drying and roasting, which is generally performed at a temperature of 135 °C for 35 to 40 minutes, are DNA destructive processing steps such that cocoa DNA cannot be used to identify and control the genetic origin of fermented beans, roasted beans and chocolate.

Consequently, an object of the present invention resides in providing methods for the identification of the genetic material of cocoa in processed cocoa beans and chocolate.

In the studies leading to the present invention it has now surprisingly been found that contrary to the common belief of the state of the art DNA identification techniques are actually suitable to detect genetic material in processed cocoa beans and chocolate and that these techniques may efficiently be used to identify the genetic origin of cocoa.

The above object has therefore been achieved by using DNA identification techniques, in particular the techniques of PCR, RAPD, RFLP or microsatellite



for the determination of genetic material of cocoa, wherein the material to be determined is fermented beans, roasted beans or chocolate.

By using the above mentioned techniques the identification of cocoa (varieties) and its (their) origin is now feasible, e.g. from which particular variety the cocoa beans contained in an end product, such as chocolate, have been derived from. This option will also represent an appealing feature so as to determine whether some particular cocoa, such as e.g. cocoa that has been modified or improved by common techniques or that has subject to genetic engineering, has been utilized in the preparation of the chocolate of interest.

Moreover, the present invention allows for an efficient control of the origin of fermented dry beans at the stage of export from/import into a particular country, for an efficient control of the origin of roasted beans in the factory during their processing and for an efficient control of chocolate on the shelf.

The techniques of PCR (Ehrlich et al., supra), RAPD (Random Amplified Polymorphic DNA; Welsh et al., "Fingerprinting genomes using PCR with arbitrary primers", *Nucleic acids Res.* **18** (1990), 7213-7218), RFLP (restriction fragment length polymorphism; Botstein D. et al., "Construction of a genetic map in man using restriction fragment length polymorphisms" *Am. J. Hum. Genet.* **32** (1980), 314-331) and microsatellite identification (Tautz, D. et al, "Hypervariability of simple sequences as a general source for polymorphic DNA markers, *Nucl. Acid. Res.* **17** (1989), 6463-6471; Weber J. at al., "abundant class of human DNA polymorphisms which can be typed using the polymerase chain reaction", *Am J. Hum. Genet.* **44** (1989), 388-396) are known and described in the art. The above mentioned documents are incorporated herewith by way of reference.

In a preferred embodiment the DNA to be detected according to the above techniques is a DNA derived from the 5 S gene of cocoa or the SSP gene (Seed Storage Protein gene; Spencer E. et Hodge, R., *Planta* **186** (1992, 567-576).

In a further preferred embodiment the DNA to be detected is derived from mitochondrial and/or chloroplastic DNA. This DNA shows the advantage that it

is stably propagated from mother to daughter trees thus enabling to assign the fruits of descendent trees to a particular original tree/variety. This option also applies when the descendents are derived from one female and different male trees.

By applying DNA identification techniques, preferably PCR, RAPD, RFLP or microsatellite identification techniques, with fermented and/or roasted cocoa beans and/or chocolate it will likewise be possible to identify the presence of a raw material derived from an variety, that has been modified or improved by common breeding techniques or that has been genetically modified, in the end product.

In the figures:

Fig. 1A shows a picture of an agarose gel stained with ethidium bromide, on which cocoa DNA from different sources have been applied.

Fig. 1B shows an autoradiography of a hybridisation experiment, wherein the DNA from the agarose gel of Fig. 1A has been transferred to a nylon membrane and has been hybridized with radioactively labelled cocoa DNA using  $^{32}\text{P}$  as the radioactive isotope.

Fig. 2 shows a picture of an agarose gel stained with ethidium bromide, in which the 5S intergenic spacer on different cocoa samples has been amplified via PCR.

Fig. 3 shows a picture of an agarose gel stained with ethidium bromide, wherein the intron 1 and the exon 2 of the SSP gene has been amplified via PCR on different samples.

Fig. 4 shows RAPD profiles from different cocoa samples.

In the following examples the techniques for the determination of cocoa DNA in fermented and/or roasted beans and chocolate will be described. It is,

however, understood that the present invention is not limited to the examples but is rather embraced by the scope of the appended claims.

Unless otherwise indicated the techniques applied are as described in Sambrook J., Fritsch E.F., Maniatis T (1989) Molecular cloning, A Laboratory Manual, second edition, Cold Spring Harbor Laboratory Press (USA).

### Example 1

DNA extraction from fermented cocoa beans:

400 mg of cocoa embryo axes were collected, ground in liquid nitrogen and homogenised in 10 ml of extraction buffer (Sorbitol 0.35 M, Tris-HCl 0.1 M pH 8, EDTA 5 mM, Sodium bisulphite 5g/l) using an IKA-Ultra-Turrax T25 (Janke & Kunkel). The resulting homogenate is filtrated through Miracloth (Calbiochem) and centrifuged (15 min, 1000 g, 4°C). The pellet is resuspended in 5 ml of extraction buffer and centrifuged (15 min, 1000 g, 4°C). The pellet is resuspended in 0.5 ml of extraction buffer with 0.5 ml of lysis buffer (Tris-HCl 0.2 M pH 8, EDTA 50 mM, NaCl 2 M, CTAB 2%) and 0.1 ml of 5% Sarkosyl and incubated for 60 min at 65°C. The lysate was then extracted with 1 ml of chloroform/isoamyl alcohol (24/1:v/v). After centrifugation (15 min at 10.000 g, 4°C) the DNA in the aqueous phase was precipitated by adding one volume of isopropanol and standing overnight at -20°C. The DNA was collected by centrifugation, resuspended in 1 ml of TE buffer (Tris-HCl 10 mM pH 8, EDTA 0.1 mM) and purified with RNase A (Sigma) at a final concentration of 200 µg/ml at 37°C for 30 min, and Proteinase K (Boehringer Mannheim) at a final concentration of 400 µg/ml at 56°C for 1 h. Proteins were removed with phenol/chloroform/isoamyl alcohol (25/24/1:v/v/v) treatment. After centrifugation (15 min, 1000g, 4°C) the DNA from the aqueous phase was precipitated overnight with a mixture of 3 M sodium acetate, pH 5.2, and cold ethanol (0.1/2:v/v; -20°C). The DNA was hooked on a glass rod, washed in 70% ethanol solution and resuspended in 50 µl of TE buffer.

### Example 2

DNA extraction from roasted cocoa beans or nibs and chocolate:

60 g of dark chocolate (Nestlé, Grands chocolats, Noir (74% cocoa); Auchan, Vendome Noir (52% cocoa), freely available on the market or 25 g of nib (from Nestlé St. Menet, France) were ground in liquid nitrogen and subsequently homogenised in 200 ml of extraction buffer as previously described. The homogenate was filtrated through Miracloth and centrifuged for 15 min at 1000 g at 4°C. The pellet was resuspended in 50 ml of extraction buffer and centrifuged (15 min, 1000 g, 4°C). The pellet was resuspended in 5 ml of extraction buffer with 7 ml of lysis buffer and 2.4 ml of 5% Sarkosyl and incubated for 60 min at 65°C. The resulting lysate was then extracted with 10 ml of chloroform/isoamyl alcohol (24/1:v/v). After centrifugation (15 min, 1000 g, 4°C) the DNA in the aqueous phase was precipitated by adding one volume of isopropanol and standing overnight at -20°C. After two successive centrifugations (15 min at 10.000 g, 4°C) the resulting pellets were pooled and resuspended in 5 ml of TE buffer. Two additional purification steps with phenol/chloroform/isoamyl alcohol (25/24/1:v/v/v) were realised and the DNA from the aqueous phase finally obtained was precipitated overnight by adding 3 M sodium acetate pH 5.2 and cold ethanol (0.1/2:v/v) and standing at -20°C. After centrifugation (15 min, 1000 g, 4°C) the DNA was resuspended in 30 µl of TE buffer. The DNA was purified with RNase A (final concentration: 200 µg/ml) at 37°C for 10 min.

### Example 3

PCR DNA amplification:

Two types of PCR amplification were tested, first on repeated DNA sequences from rDNA 5S gene (1000 to 50000 copies per haploid genome) and second on a single gene coding for a cocoa seed storage protein (SSP) with 3 to 5 copies per haploid genome.

DNA amplification for the intergenic spacer of ribosomal 5S gene:

Purified DNA obtained from cocoa beans, nib or chocolate as illustrated above was diluted in TE buffer with a ratio of 1/500. The DNA amplification is performed with 5  $\mu$ l of each diluted DNA samples and 45  $\mu$ l of the PCR buffer containing 200  $\mu$ M of each dNTP, 25 ng of each primer and 2 U of *Taq* polymerase (Stratagène) with a final concentration of 1.5 mM magnesium chloride.

The DNA sequences of the two primers used are:

5'-TTTAGTGCTGGTATGATCGC-3' (SEQ. ID. No. I)

5'-TGGGAAGTCCTCGTGTTGCA-3' (SEQ. ID. No. II)

These two primers were designed according to Kolchinsky et al., "Portraying of plant genomes using polymerase chain reaction amplification of ribosomal 5S genes", Genome **34** (1991), 1028-1031. Amplifications were performed in a Bio-Med Thermocycler 60/2 (B. Braun) programmed for a preliminary 2-min denaturation at 94°C and 30 cycles including 1 min at 94°C, 1 min at 55°C and 1 min at 72°C. A final DNA extension cycle was performed at 72°C for 7 min.

DNA amplification for the seed storage protein (SSP) gene:

The DNA amplification procedure was the same as previously described with the proviso that the two DNA primers designed to detect the presence of the SSP gene are as follows:

5'-GGCAATTTACTTCGTGACAAACG-3' (SEQ. ID. No. III)

5'-CTCATATTTGCCAGGAGAATTAAC-3' (SEQ. ID. No. IV)

Example 4

Random amplified polymorphic DNA (RAPD):

DNA amplification was performed as already described above with the proviso that the primers used are decamers randomly designed from Operon Technologies, Alameda, California.

The Codes for the corresponding primers are as follows:

AG15	5'-CCCACACGCA-3'	(SEQ ID. NO V)
AM10	5'-CAGACCGACC-3'	(SEQ ID NO. VI)
Z06	(to be found in the PCR kit from Operon)	

Amplifications were programmed for 45 cycles including 1 min at 94°C, 1 min at 37°C and 2 min at 72°C. A final DNA extension cycle was performed at 72°C for 7 min according to Williams et al., "DNA polymorphisms amplified by arbitrary primers are useful as genetic markers", Nucl. Acids Res. **18** (1990), 6531-6535.

DNA electrophoresis and hybridisation:

The different DNA fragments were separated on a 1.4 % agarose gel followed by staining with ethidium bromide. Southern blots were made according Southern (Sambrook, supra) using alkaline capillary transfer to nylon membranes (Appligène). The DNA probes were labelled using <sup>32</sup>P dCTP with the Megaprime kit (Amersham) and were hybridised on the membranes overnight at 65°C in hybridisation buffer (5% w/v SDS, 5 X SSC, 5 X Denhart's solution, 40µg/ml of heterologous DNA). Post-hybridisation treatments consisted of three high- stringency washes (2 X SSC, 0.1% SDS; 1 X SSC, 0.1% SDS; 0.2 X SSC, 0.1% SDS) each at 65°C for 30 min.

### Example 5

DNA detection from cocoa bean to chocolate:

The protocol used allowed the detection of DNA of high molecular weight in fermented cocoa beans and chocolate. As negative control coffee and hazelnut DNA were used. To ensure that the DNA extracted from nib and chocolate was derived from cocoa the DNA transferred to the nylon membrane has been hybridized with radioactively labelled cocoa total DNA purified from cocoa tree leaves. The autoradiography (Figure 1B) indicates positive homology for the cocoa samples from leaves to chocolate which demonstrates that the DNA purified from these samples was originated from cocoa.

PCR DNA amplification from cocoa samples:

In order to test the possibility to detect specific DNA sequences from the different samples specific DNA primers were used to amplify first the intergenic spacer from the 5S ribosomal gene tandemly repeated and second a part of the seed storage protein gene (SSP).

The 5S intergenic spacer amplifications were successful on all the DNA samples tested from leaves to chocolate (Figure 2). These PCR DNA amplifications were detected as a faint DNA smear indicating a substantial degradation of the DNA during the processing steps of the cocoa beans with discrete bands to be detected from approximately 160 bp to up than more 1000 bp for the leaf and fresh bean samples. These multiple DNA bands amplification were the result of the tandem repetition organisation of the target 5S ribosomal gene.

The part of the SSP gene amplified corresponds to the intron 1 and exon 2 thereof and yields a PCR product of 312 bp in length which may be detected in all the samples tested (Figure 3). This result allows to consider that it is possible to detect a single DNA gene copy from leaf to chocolate since SSP gene represents a low gene copy number in *Theobroma cocoa* haploid genome (3 to 5 copies).

DNA polymorphism was detected for the PCR DNA amplification from part of SSP gene.

#### RAPD tests on cocoa samples:

The goal of this experiment was to demonstrate that it could be possible to get RAPD fingerprint from processed cocoa samples to determine the genetic origin(s) of the raw material used for a chocolate. Three primers were selected from Operon kits and the results are illustrated by the figure 4. DNA amplifications were obtained in two different chocolate blends (Nestlé, Grands chocolats, Noir (74% cocoa); Auchan, Vendome Noir (52% cocoa)) showing that it is indeed possible to relate the end product to particular raw materials.



### Claims


1. Use of a DNA detection technique for the determination of genetic material of cocoa in fermented and/or roasted beans and/or chocolate.
2. The use according to claim 1, wherein the DNA detection technique is selected from the group consisting of PCR, RAPD, RFLP or microsatellite identification.
3. The use according to any of the claims 1 or 2, wherein the DNA to be detected is derived from the 5 S gene of cocoa.
4. The use according to claim 1, wherein the DNA to be detected is derived from the SSP gene of cocoa.
5. The use according to claim 1, wherein the DNA to be detected is derived from mitochondrial and/or chloroplastic DNA of cocoa.
6. The use according to any of the preceding claims, wherein the cocoa from which the fermented and/or roasted beans and the chocolate are derived, is a variety, that has been modified by common breeding techniques or that has been modified by genetic engineering.

THIS PAGE BLANK (USP 10)

### Summary

The present invention pertains to the use of DNA detection techniques for the determination of cacao in fermented and/or roasted beans and chocolate varieties and/or cacao varieties that have been modified by common breeding techniques or that have been modified by genetic engineering. In particular the DNA detection techniques are selected from the group consisting of PCR, RAPD, RFLP or microsatellite identification.

THIS PAGE BLANK (USPTO)

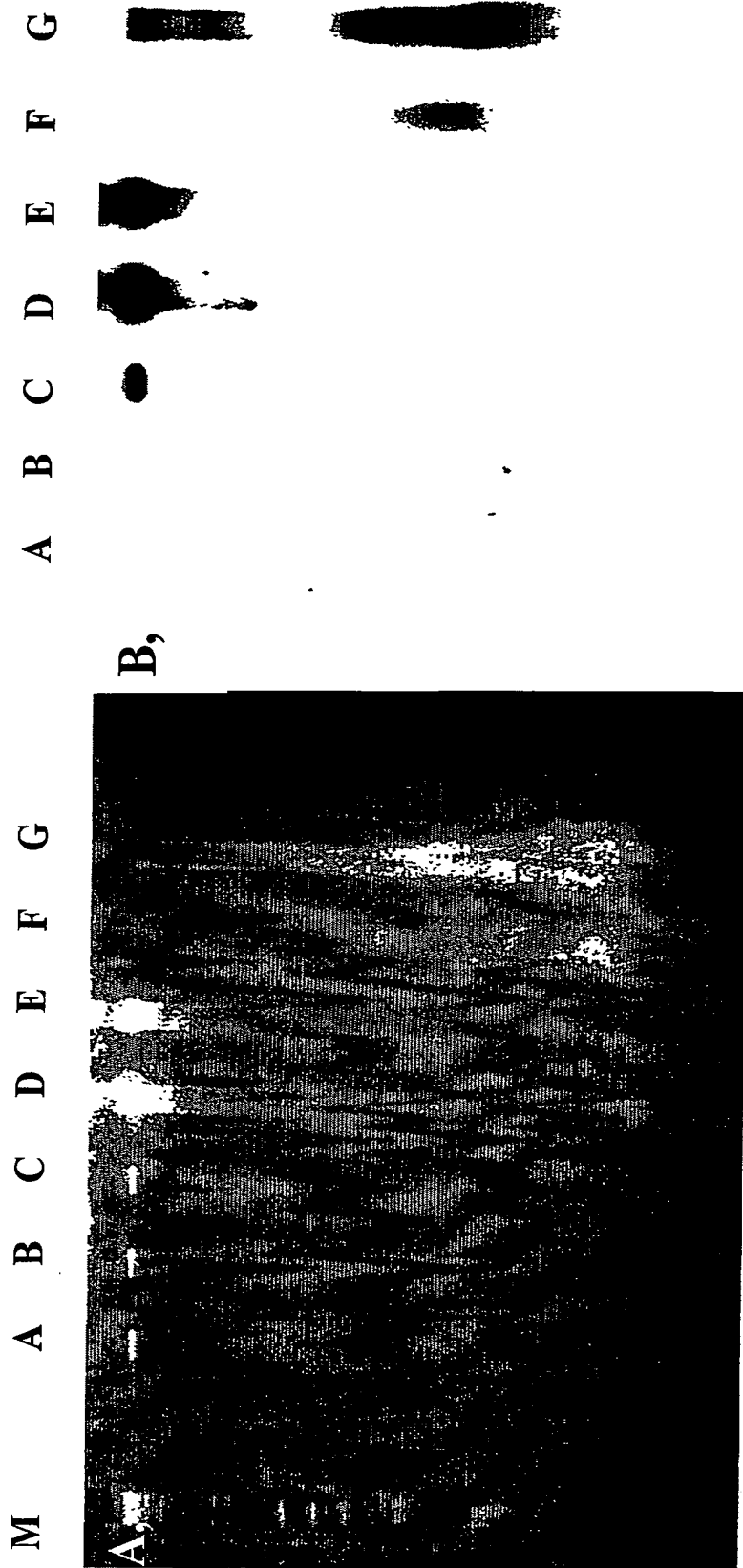


**Figure 1. DNA detection in various cacao and cocoa samples**

1A, Native DNA detection in agarose gel using ethidium bromide


1B, Hybridisation of total cacao DNA on membrane transfer of native DNA agarose gel

M: indicates molecular size marker ( $\lambda$ HindIII and  $\Phi$ 174/HaeII), A is a DNA control from coffee leave, B is a DNA control from hazelnut leave, C is a DNA control from cacao leave, D is DNA sample from fresh cacao seed embryo, E is a DNA sample obtained with fermented cacao beans, F is a DNA sample from roasted nib and G is a DNA sample from dark chocolate (Nestlé Noir).

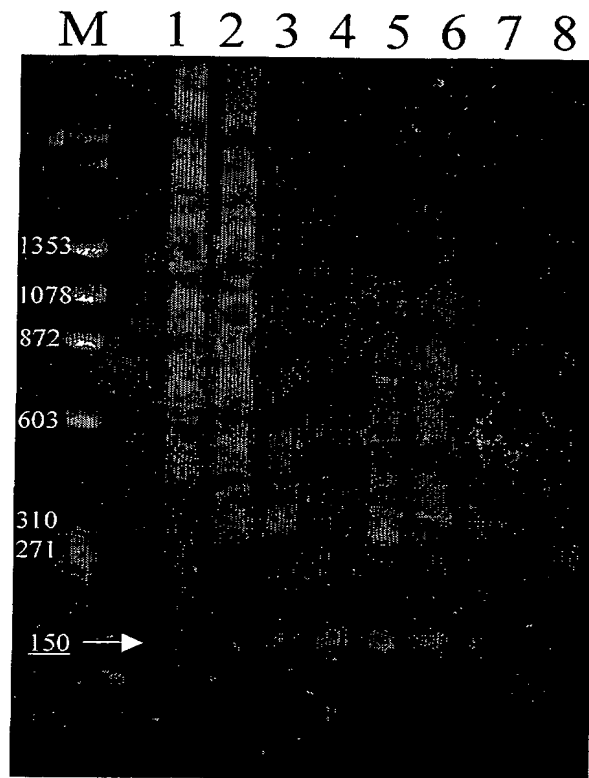


**Figure 2: PCR-DNA amplification of 5S intergenic spacer on different cacao samples.**

M: indicates molecular size marker in base pairs ( $\lambda$ HindIII and  $\Phi$ 174/HaeIII), 1: Cacao leaves, 2: Cacao fresh bean, 3: Cacao fermented bean, 4 & 5: Cocoa roasted nib, 6 & 7 dark chocolate (Nestlé Noir), 8: negative control

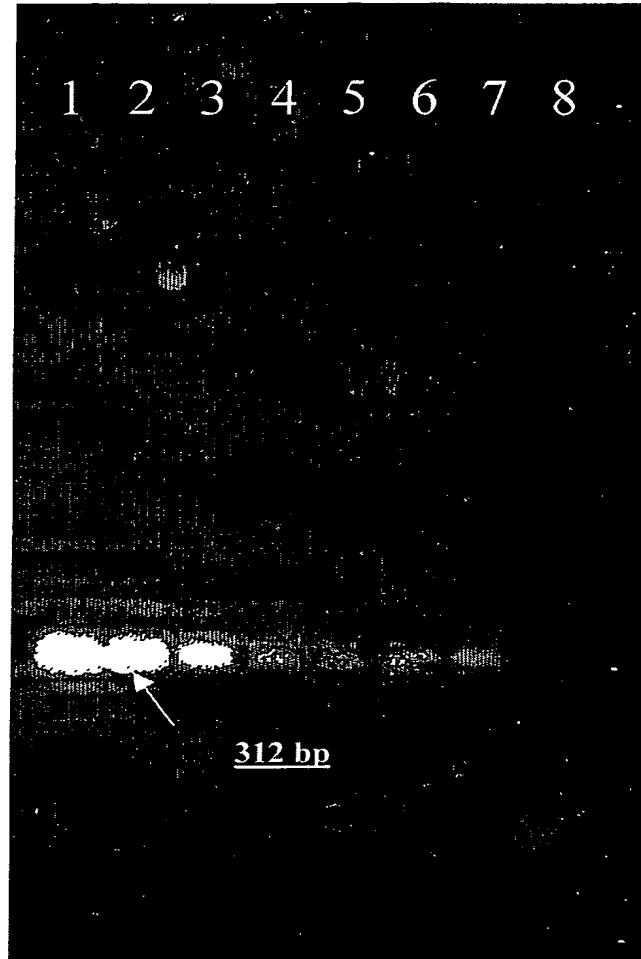

**Figure 3: PCR-DNA amplification of intron 1 and exon 2 of Seed Storage Protein gene (SSP)**  
1: Cacao leaves, 2: Cacao fresh bean, 3: Cacao fermented bean, 4 & 5: Cocoa roasted nib, 6 & 7 dark chocolate (Nestlé Noir), 8: negative control

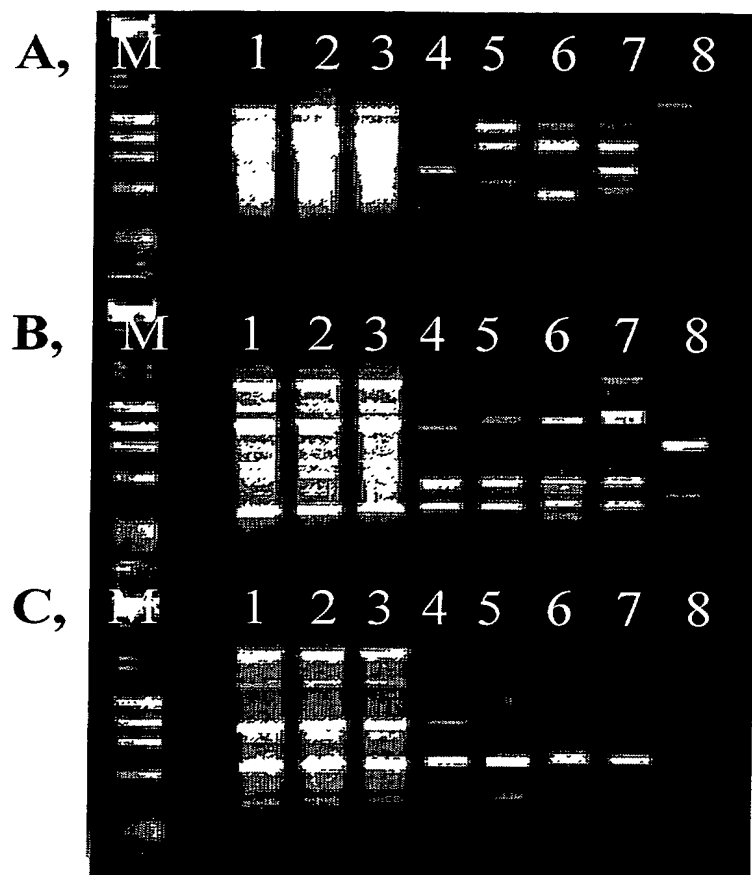

**Figure 4: RAPD profiles from various cacao and cocoa samples**

A: Z06 primer, B: AG15 primer, C: AM10. M: indicates molecular size marker ( $\lambda$ HindIII and  $\Phi$ 174/HaeIII), 1, 2 and 3 are cacao leaf samples, 4 and 5 are cocoa samples from "Nestlé Noir", 6 and 7 are cocoa from "Vendome" and 8 indicates the negative control





**Nestlé®**  
*Grands Chocolats*

# NOIR

CHOCOLAT NOIR DE DÉGUSTATION A 74% DE CACAO

CHOCOLAT NOIR



A conserver dans un lieu sec et frais  
(10° à 18° C) hors du réfrigérateur.

**100 g e**

A CONSERVER DANS UN LIEU SEC ET FRAIS, A L'ÉLOIGNER DE LA CHALEUR.

**Nestlé® Grands Chocolats NOIR**

**Nestlé® Grands Chocolats**  
et  
**LES PLUS GRANDES ÉMOTIONS DE LA MUSIQUE**

Nestlé Grands Chocolats a créé pour vous une collection exclusive de 3 CD (à titre) :  
Musique Classique et des compositeurs illustres : B.A.B., Chopin, Mozart...  
Variété Française et des chansons françaises (enregistrement original) : "La Ballade des gens  
bonnes", "Qu'on se papillonne", "Félicité"...  
Jazz et les plus grands du jazz américain : Louis Armstrong, Billie Holiday, Glenn Miller...  
Pour recevoir le CD de votre choix, envoyez le bon à Nestlé Grands Chocolats :  
1 CD à 9 francs 90 (hors taxes) et 10 francs 90 (hors taxes) en plus.

« Collectionnez 4 codes-barres à découper au dos de chacune des tablettes Nestlé Grands Chocolats  
trouvées dans les tablettes fourrées Nestlé Grands Chocolats ou sur papier libre.  
• Joignez le ou les tickets de caisse où figurent vos achats de tablettes avec 8 timbres au  
taux "lettre" en vigueur.  
• Collez vos codes-barres et inscrivez le CD de votre choix : Musique Classique ou Variété  
Française ou Jazz, avec vos coordonnées complètes sur la cartonne collector que vous  
trouverez dans les tablettes fourrées Nestlé Grands Chocolats ou sur papier libre.

CADRE : Votre 1er code-barre vous est déjà offert sur la cartonne collector à l'adresse des  
collecteurs fourrés Nestlé Grands Chocolats (rte d'Amboise, Trélon, Druis de Mennecy, Orléans).  
Envoyez le tout à : "Les Plus Grandes Émotions de la Musique", Code 5037 - 95400 Paris Clichy  
avant le 31.12.1994. Vous recevrez le CD de votre choix dans un délai de 4 semaines.  
\*Offre valable jusqu'au 31.12.1994. En cas de rupture de stock, Nestlé pourra se réserver la possibilité d'accepter un autre CD.

Nestlé Grands Chocolats  
VOUS OFFRE

LES PLUS  
GRANDES  
ÉMOTIONS  
DE LA  
MUSIQUE

Recevez votre CD  
des meilleurs titres de  
musique classique,  
de variété française  
ou de jazz


Pour de plus amples  
informations sur nos chocolats  
SERVICE CONSOMMATEURS  
NESTLÉ FRANCE  
BP 920 NOISY  
77146 MARNE-LA-VALLÉE Cedex 2 FRANCE



Nestlé Grands Chocolats  
bénéficie du savoir-faire  
et du souci de qualité de Nestlé,  
pour des chocolats  
plus forts en plaisir

Pour avoir  
cette partie

Annex 1: "Nestlé Noir" blend used to detect DNA in dark chocolate




CHOCOLAT

**Vendôme**

**NOIR**

**Vendôme**



Le chocolat noir VENDÔME a été élaboré à partir des meilleurs fèves des Tropiques. Fabriqué dans le respect de la tradition chocolatière, c'est un mariage réussi entre finesse du chocolat, force et équilibre du cacao. D'une saveur raffinée, le chocolat VENDÔME est un plaisir de haute gourmandise qui régalera gourmets et gourmands.

Valores energéticas y nutricionales por 100 g. / Valeurs énergétiques et nutritionnelles par 100 g. / Valores energéticos y nutricionales por 100 g.

<b>2260 kJ</b> (540 kcal)	<b>Protéinos</b> Proteínas Proteínas <b>5,7g</b>	<b>Glúcidos</b> Maltosa de carbono Glucosa Glucosa <b>58,5g</b>	<b>Lípidos</b> Grasas Lípidos Lípidos <b>31,7g</b>
------------------------------	---	--	--


**Chocolat noir.**  
INGREDIENTS: Liqueur de cacao, sucre, beurre de cacao, lactose, arôme naturel, lecithine de soja, huile de palme, arôme. Cacao 52 % minimum.  
Conservation: à l'abri de l'humidité. Température idéale de conservation: 18-20°C.

**Chocolato negro.**  
INGREDIENTS: Liqueur de cacao, azúcar, manteca de cacao, lactosa, aroma natural, aceite de soja, aceite de palma, aroma. Cacao 52 % mínimo.  
Almacenaje: en un lugar seco y fresco. Temperatura ideal de conservación: 18-20°C.

**Chocolato fondente.**  
INGREDIENTS: Liqueur de cacao, sucre, beurre de cacao, lactose, émulsifiant: lecithine de soja, huile de palme, arôme. Cacao 52 % minimum.  
Conservation: à l'abri de l'humidité. Température idéale de conservation: 18-20°C.

**Chocolato preto.**  
INGREDIENTS: Liqueur de cacao, açúcar, manteiga de cacao, lactose, emulsificante: lecitina de soja, óleo de palma, aroma. Cacao 52 % mínimo.  
Armazenagem: em local seco e fresco. Temperatura ideal de conservação: 18-20°C.

Made in U.S.A.  
Poids net (Peso neto) / Peso neto (Peso líquido): **200g e**



3 254560 888098

CHOCOLAT noir. Appréhender le chocolat noir dans l'absolu, c'est un plaisir. C'est la haute cuisine du chocolat. Vendôme, c'est la tradition. Appréhender le chocolat noir dans l'absolu, c'est un plaisir. C'est la haute cuisine du chocolat. Vendôme, c'est la tradition.

**SANTÉFANT**  
DÉLICIEUSE  
CHOCOLAT  
SANTÉFANT PERSONALIZZATO  
GUSTO AROMATICO - 85-85-85  
100% MELLINONE D'ASCO CENNI

Annex 2: "Vendome" blend used to detect DNA in dark chocolate

5. Nov. 1998

Sequence Listings

## (1) GENERAL INFORMATION:

- (i) APPLICANT: Societe des Produits Nestlé, S.A.  
P.O. Box 353  
CH-1800 Vevey  
INVENTORS: Vincent Petiard  
David Crouzillat

TITEL OF THE INVENTION: Use of DNA identification techniques  
for the determination of genetic material  
of cocoa in fermented or roasted beans  
and chocolate

(iii) NUMBER OF SEQUENCES : 6

## (iv) ADDRESS:

- (A) Address: Kirschner & Kurig  
(B) Street: Sollnerstr. 38  
(C) City: Munich  
(D) State: Bavaria  
(E) Country: Germany  
(F) Regional Code: 81479

## (v) COMPUTER READABLE FORM:

- (A) Medium type: Floppy Disk  
(B) Computer: IBM PC Compatibel  
(C) System: PC-DOS/MS-DOS  
(D) Software: Patentin Release Nr. 1, Version Nr. 1,25

## (vi) INFORMATION OF THE APPLICATION:

- (A) Application number: to be obtained  
(B) Application date: herewith  
(C) Classification:

(vii) PRIORITY DATES: none

## (viii) Attorney/Agent Information

- (A) Name: Dr. Straus, Alexander  
(B) Registration number: 85880  
(C) Reference : E 80,022 EP

## (ix) TELECOMMUNICATION:

- (A) Telephone: (089) 749 858 - 0  
(B) Telefax: (089) 749 858 - 11

## (2) INFORMATION FOR SEQ ID. NO. 1:

## (i) SEQUENCE CHARACTERISTICS:

- (A) length: 20 nucleotides
- (B) type: nucleic acid
- (C) strandedness: single stranded
- (D) topology: linear

## (ii) MOLECULE TYPE: DNA

## (xi) SEQUENCE DESCRIPTION: SEQ ID. NO.1:

TTTAGTGCTG GTATGATCGC

20

## (3) INFORMATION FOR SEQ ID. NO. 2:

## (i) SEQUENCE CHARACTERISTICS:

- (A) length: 20 nucleotides
- (B) type: nucleic acid
- (C) strandedness: single stranded
- (D) topology: linear

## (ii) MOLECULE TYPE: DNA

## (xi) SEQUENCE DESCRIPTION: SEQ ID. NO.2:

TGGGAAGTCC TCGTGTTGCA

20

## (4) INFORMATION FOR SEQ ID. NO. 3:

## (i) SEQUENCE CHARACTERISTICS:

- (A) length: 23 nucleotides
- (B) type: nucleic acid
- (C) strandedness: single stranded
- (D) topology: linear

## (ii) MOLECULE TYPE: DNA

## (xi) SEQUENCE DESCRIPTION: SEQ ID. NO.3:

## (5) INFORMATION FOR SEQ ID. NO. 4:

## (i) SEQUENCE CHARACTERISTICS:

- (A) length: 20 nucleotides
- (B) type: nucleic acid
- (C) strandedness: single stranded
- (D) topology: linear

(ii) MOLECULE TYPE: DNA

(xi) SEQUENCE DESCRIPTION: SEQ ID. NO.4:

CTCATATTTG CCAGGAGAAT TAAC

24

(6) INFORMATION FOR SEQ ID. NO. 5:

(i) SEQUENCE CHARACTERISTICS:

(A) length: 10 nucleotides

(B) type: nucleic acid

(C) strandedness: single stranded

(D) topology: linear

(ii) MOLECULE TYPE: DNA

(xi) SEQUENCE DESCRIPTION: SEQ ID. NO.5:

CCCACACGCA

10

(7) INFORMATION FOR SEQ ID. NO. 6:

(i) SEQUENCE CHARACTERISTICS:

(A) length: 10 nucleotides

(B) type: nucleic acid

(C) strandedness: single stranded

(D) topology: linear

(ii) MOLECULE TYPE: DNA

(xi) SEQUENCE DESCRIPTION: SEQ ID. NO.6:

CAGACCGACC

10

**THIS PAGE BLANK (USPTO)**

**This Page is Inserted by IFW Indexing and Scanning  
Operations and is not part of the Official Record**

**BEST AVAILABLE IMAGES**

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- ☐ BLACK BORDERS
- ☐ IMAGE CUT OFF AT TOP, BOTTOM OR SIDES
- ☐ FADED TEXT OR DRAWING
- ☐ BLURRED OR ILLEGIBLE TEXT OR DRAWING
- ☐ SKEWED/SLANTED IMAGES
- ☐ COLOR OR BLACK AND WHITE PHOTOGRAPHS
- ☐ GRAY SCALE DOCUMENTS
- ☒ LINES OR MARKS ON ORIGINAL DOCUMENT
- ☐ REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY
- ☐ OTHER: \_\_\_\_\_

**IMAGES ARE BEST AVAILABLE COPY.**

**As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.**

THIS PAGE BLANK (USPTO)